


成品尺寸 48x148mm

Multifunctional Eye Beauty Device

USER MANUAL

Eye Beauty Device

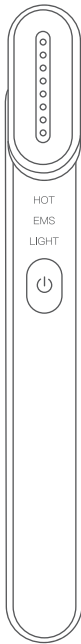

Before using this product, please read the "User Manual" carefully for reference and keep it properly.

Note: If the model in this manual may differ from the actual product due to the iteration of the product appearance, please refer to the actual product.

Product Introduction

Combined with LED photon rejuvenation, 38-42 °C hot compress, EMS micro current stimulation, this product with about 12000 times/minute of micro vibration is used after applying cream (or essence) and other skin care products, which can promote the absorption of skin care products nutritional essence, help reduce dark circles/bag /puffiness/crow eyes, fade wrinkles, and delay the formation of true wrinkles.

Red Light

Red light can improve cell activity, accelerate metabolism and promote collagen formation, which can effectively whiten and lighten spots, rejuvenate skin and remove wrinkles, increase skin elasticity, improve skin condition, repair damaged skin, minimize pores, tighten skin, anti-aging, antioxidant and so on.

Blue Light

Blue light has anti-bacterial and anti-inflammatory qualities, which not only effectively removes acne and sores, sterilizes and reduces inflammation, but also repairs the skin without leaving scars or pigmentation

Green Light

Green light can reduce the secretion of skin oil, balance the skin, effectively relieve stress, dredge lymph and eliminate face edema.

Yellow Light

Yellow light smoothes skin and reduces redness, fade dark spots, freckles, and blemishes.

Purple Light

Purple light is a dual-frequency light with red light and blue light so that it combines the effects of both lights, especially suitable for treating acne and acne scars.

Cyan Light

Cyan light can gradually increase cellular energy and has a good effect on metabolism.

White Light

White light can penetrate into the deeper layer of the skin, accelerating the skin's metabolism, effectively breaking down pigmentation, improving fine lines and sagging skin.

Major Function

1. LED photorejuvenation
2. High-frequency vibration to promote the absorption of essence
3. Soothing hot compress to promote the absorption of essence
4. EMS microcurrent for skin tightening

Usage Method

Operation Method

01

Apply proper amount of eye cream or essence around the eyes to keep the skin moist.



02

Slide the massage head gently, lift and massage the skin around the eyes, repeat 15-20 times;



03

Press the massage head gently to improve the absorption of eye cream or essence. Stay at each point for about 10 seconds and repeat for 2-3 times.



04

You can also massage the cheeks, chin and other parts according to the above operation methods. You can use the lifting method to make your skin better.



05

After use, wipe the massage head with a clean cloth, alcohol cotton

06

Besides eye massage, you can also massage the whole face and other parts of the body (such as the waist and legs)

Matters Need Attention

1. Insert the USB terminal into the customer's 5V USB charging head and align the Type-C charging head with the charging post on the bottom of the product, which will automatically charge regardless of direction;
2. When the indicator light fluctuates and blinks, it means it is in charging status;
3. When the indicator light is always on, it means that it is fully charged.
4. If the indicator light flashes frequently when the product is in use, it means the power is insufficient and should be charged immediately, otherwise it will stop working;
5. The charging cable should be properly stored and should not be lost.
6. The product will automatically turn off after three minutes of operation.

Warning

This product is not waterproof. Do not place the product near a bathtub sink or washbasin to prevent damage to the product due to water ingress. It should be ensured that the battery is fully charged at least once within 3 months. Otherwise, the battery may dry out and become damaged;

After-sale Service

If there is anything wrong with our product or you have any questions, please contact our customer service staff by email without hesitation.

Suggestions for Use:

- Please fully charge the battery when using it for the first time. Before use, please clean your skin, take off your makeup and watches, earrings, glasses, contact lenses and other metal accessories.
- No using during physiological periods. Do not use if you are under 18 years old or you are pregnant.
- If you have implanted medical equipment or you have received plastic surgery, please do not use this product.
- Apply essence on the care part to keep your skin moist.
- Don't stay in the same place for a long time to avoid scalding.
- If the skin is damaged, please stop using.
- After the care, you can do daily skin care and moisturizing facial mask. The equipment is not waterproof, so it is advisable for you to clean the residue on the equipment with dry cloth or paper towel to keep it dry.


Product Parameters

Product name:	Eve Beauty Device
Battery capacity:	350mAh
Product size:	6"x 0.63" x 0.43"
Light Modes:	7 kinds of light
Charging input:	5V/1A
Product weight:	35g
Main materials:	Metal +PC plastic
Charging time:	About 2 hours
Power:	1.5W

What's Included

- 1 Eye Beauty Device
- 1 USB Type C Charging Cable
- 1 User Manual
- 1 Storage Bag

Part Description

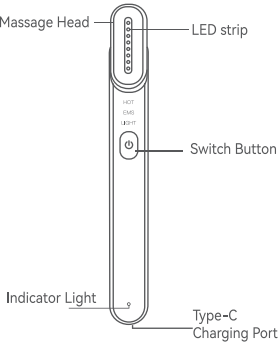

FAQ

Q: Why do I sometimes not feel the microcurrent when it touches my skin?

A: Because your skin is too dry, the machine does not sense it. You need to keep your skin moist and make the machine touch closely to your skin. It is recommended to use it together with gel or essence to achieve better results.

Q: Why does the machine sometimes fail to start?

A: It may be that the power of the machine is insufficient, it is recommended to use it after full charge.

Q: How often are they used?

A: It is recommended to use it for 5-10 minutes each time, 3-5 times a week.

Q: Why can't I use it when charging?

A: For safety. It can't be used during charging. It can only be used after charge.

Q: Is the product waterproof?

A: This product is not waterproof so that it can't be cleaned by water. It is recommended to wipe the surface of the machine with dry cloth or paper towel.

Q: How long does it take to charge?

A: It takes about 2 hours to charge, and each full charge can be used 4-5 times.